


HOW TO USE AN ASTHMA INHALER WITH A SPACER

- Carefully examine the spacer for damage, missing parts or foreign objects.




- Check the medication's expiration date.

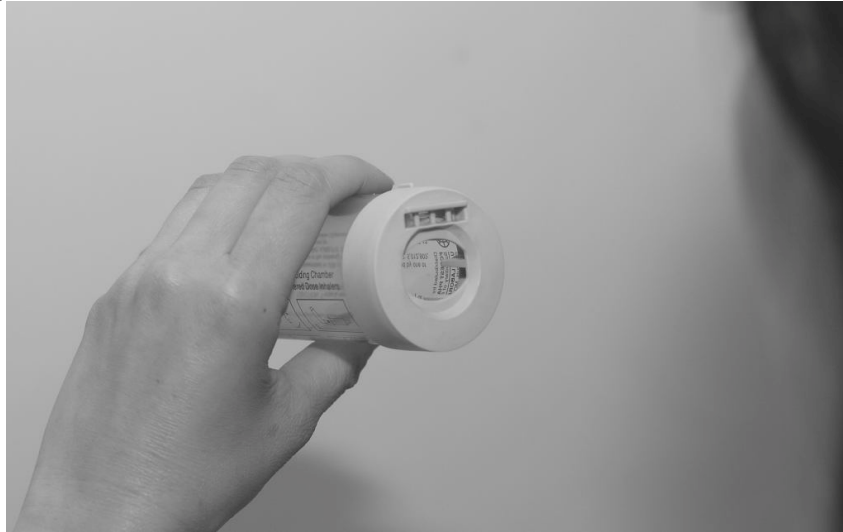

- Shake the medication's bottle well before each use



- Place the bottle in the plastic container and insert the inhaler into the back piece of the spacer.

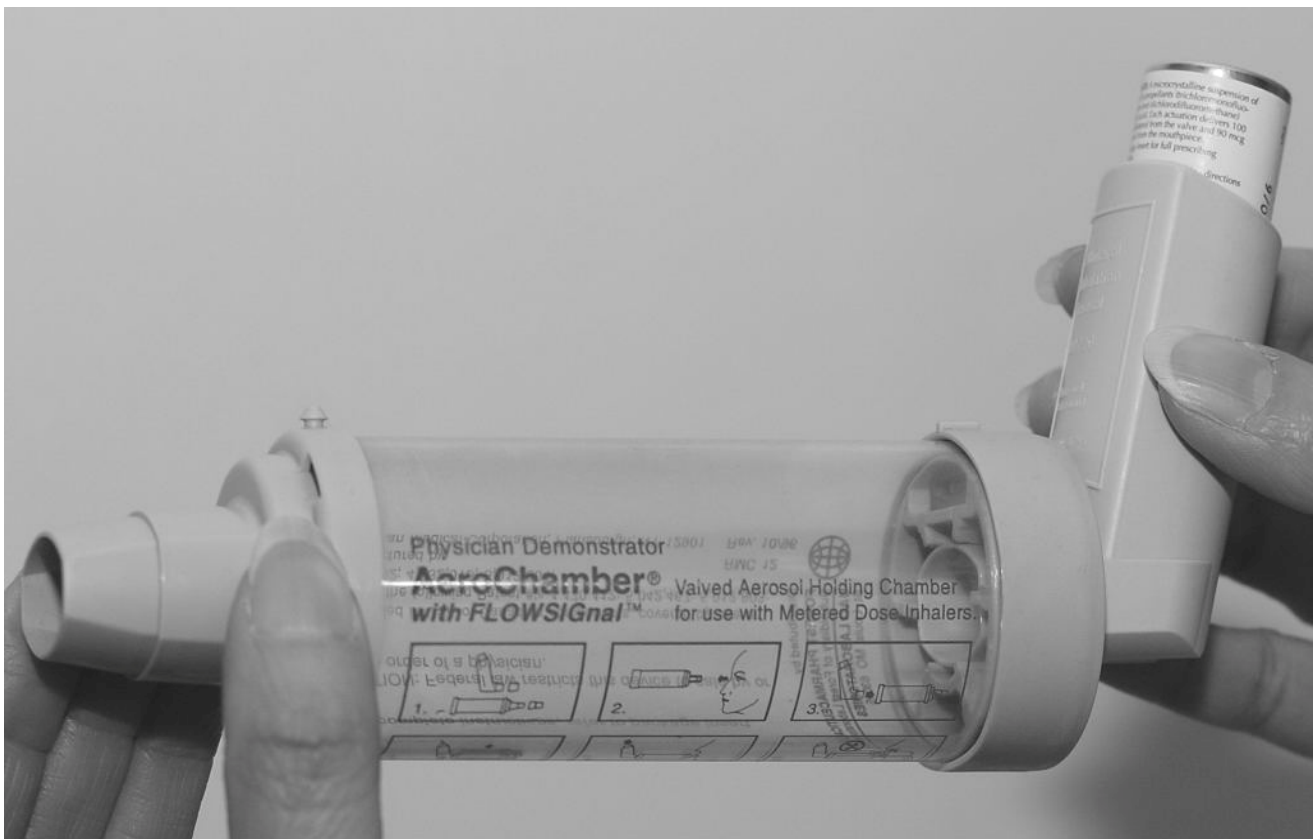

- Apply mask to face and ensure that there is a good seal. Child should stand up or be standing.



- Depress inhaler at beginning of slow inhalation.



Maintain seal for six (6) breaths after depressing inhaler. Administer one puff at a time. Do not administer more than two puffs.

Spacer Cleaning Instructions

- Clean spacer after first use, then weekly.
- Remove back piece only. Do not remove mask or tamper with valves during cleaning.



- Soak both parts for 15 minutes in warm water with liquid detergent. Agitate gently.
- Rinse in clean water.
- Shake out excess water. **Do not rub dry.**
- Let air dry in vertical position.



- Replace back piece when unit is **completely dry** and ready for use. **Do not leave the spacer unattended with children. This is a medical device, not a toy.**